


IMPORTANT INSTRUCTIONS ON HOW TO USE X-TOP CORRECTLY

1

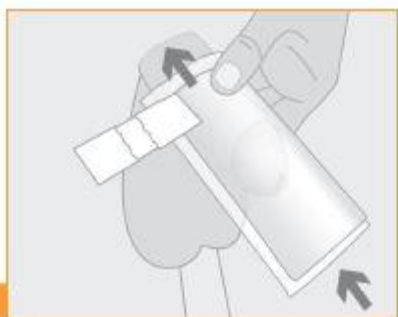

Pull the x-top pouch over the penis as far as possible.

2



Tighten and affix the hook and loop closure. Close it tight so it does not slip. The closure is adjustable so it can be re-fastened as often as needed. It's important that it stays in place.

3



Wear tight fitting underwear or briefs, no boxers please.

Change the x-top when needed, but at least once a day.

4



Dispose of used x-tops with regular household waste.

Important!

To make sure the x-top stays in place, regularly check and adjust when necessary.



Reminder!

Always wear close underwear, no boxer shorts!