


Technical sheet

PERISIZE 4

vaginal probe

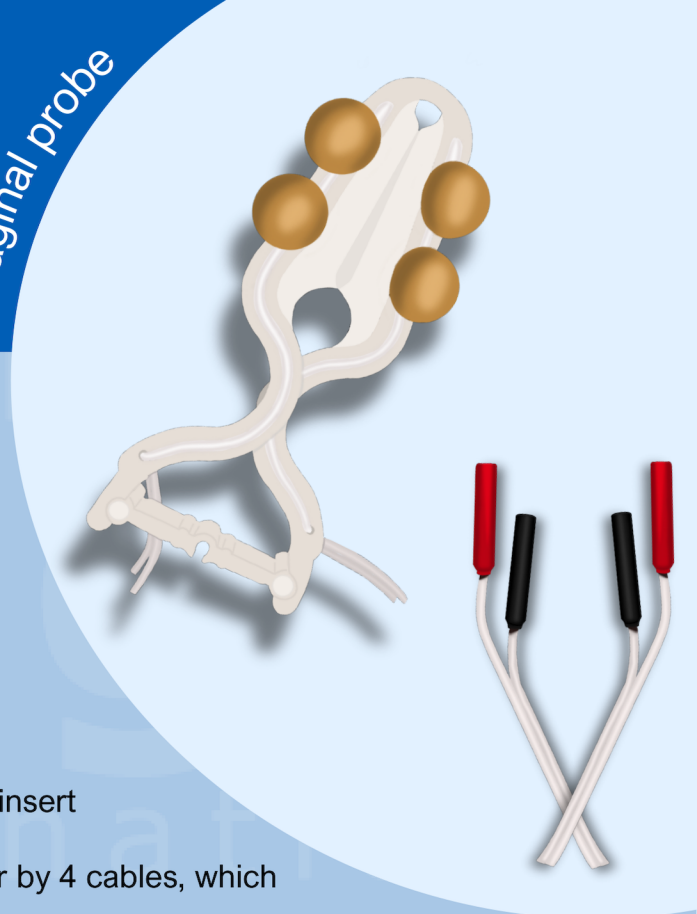

SPECIFICATIONS

- **3 hemispherical electrodes, 4 independent plugs**
- **Connection:** 2mm female banana plug
- **Diameter:** 4.8 cm
- **Length:** 11 cm
- **Weight:** 18g
- **Nickel-free**

PRESENTATION

Its ability to reduce its diameter by pinching makes it easy to insert allowing it to adapt automatically when released.

The 4 electrodes are independent and linked to the generator by 4 cables, which can be connected 2 by 2.

PRACTICAL BENEFITS OF THE PROBE

1 - Electrical stimulation - SEF

- Contemporary bilateral stimulation R & L (levator ani muscles)
- Individualized unilateral stimulation R or L (levator ani muscles)
- Bilateral stimulation front or back (pelvic floor muscles)

2 - Biofeedback - BFB/EMG

- Recording of supine or orthostatic electrical muscle activity
- Recording of R and/or L lateral electrical muscle activity in supine or orthostatic position
- Recording of individualized orthostatic stat. or dyn. global or lateral muscle activity.

THERAPEUTIC APPLICATIONS OF THE PROBE

1 - Electrical stimulation - SEF

- Direct perineal muscle deficiency
- Phasic and/or tonic perineal muscle strengthening
- Unilateral perineal analgesic treatment

- Punctual unilateral SEF if asymmetry
- Sphincter urethral insufficiency (reflex)
- Reduction of bladder hyperactivity (reflex)
- Global or lateral perineal proprioceptivity

2 - Biofeedback - BFB/EMG

- Global or lateral perineal awareness
- Detection and measurement of bilateral or unilateral muscular insufficiency

- Checking progress of bilateral or asymmetric recovery
- Control of reflex and voluntary locking during exercise, lying/standing

